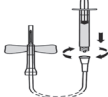
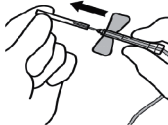
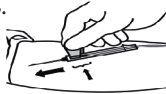
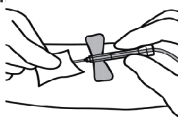
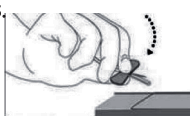


NEEDLE INSTRUCTION SHEET

General Instructions:

- Should be performed wearing gloves.
- Prime set is necessary.
- Make sure blood is present in the chamber before administration of medication.
- Dispose of the Luer cap, needle and all materials used in appropriate disposal container.

PROCEDURE for Short Term IV Administration

1. 	▪ Remove Luer adaptor with attached needle then attach syringe to the Female Luer Adaptor.									
2. 	▪ Remove needle sheath and prime set in accordance with recommended procedure. Prefill the system with medication and ensure that there is no air in the system. <table border="1"><thead><tr><th colspan="3">Priming Volume</th></tr><tr><th>Gauge</th><th>Tubing (mm)</th><th>Volume (mL)</th></tr></thead><tbody><tr><td>25</td><td>230</td><td>0.34</td></tr></tbody></table>	Priming Volume			Gauge	Tubing (mm)	Volume (mL)	25	230	0.34
Priming Volume										
Gauge	Tubing (mm)	Volume (mL)								
25	230	0.34								
3. 	▪ Insert needle in vein by holding wings. Administer over a period of 15 to 30 seconds.									
4. 	▪ Withdrawing Surshield safety winged infusion set: Apply pressure to site using your fingers as shown. ○ Grasp the translucent orange safety shield grip area with the thumb and index finger with one hand and with opposite hand, grasp tubing between thumb and index finger. ○ Cover the puncture site with a sterile gauze pad and apply pressure.									
5. 	▪ Surshield safety winged infusion set Shield Activation: One-hand Technique ○ With the wing and shield between your thumb and index finger, press the safety shield against a hard surface such as a bedside table until an audible click is heard.									
6. 	▪ Dispose of used infusion set promptly in appropriately labeled container to prevent reuse or accidental injury.									

GENERAL PRECAUTIONS:

- Handle all biologic samples and blood collection "sharps" (lancets, needles, luer adapters, and blood collection sets) in accordance with the policies and procedures of your facility.
- Obtain appropriate medical attention in the event of any exposure to biologic samples (e.g. through a puncture injury) since samples may transmit viral hepatitis, HIV (AIDS) or other infectious diseases.
- Utilize any built-in used needle protector if the blood collection device provides one.
- Discard all blood collection "sharps" in appropriate biohazard disposal containers.